

Microdispensing Instruments


Grow with your application

From compact entry level R&D applications to high-speed manufacturing




Same software for all spotters

Easily control your instruments with the intuitive and user-friendly InDot software. No programming language - just drag and drop tiles to form a sequence and loops.

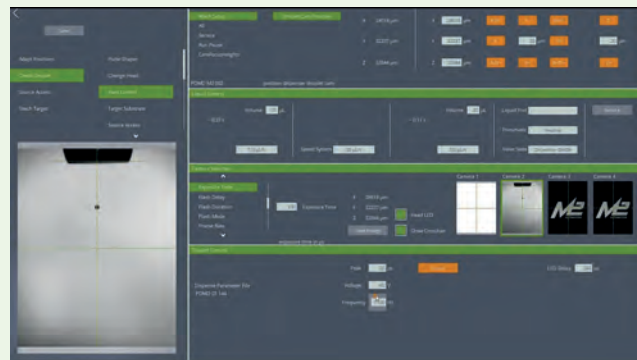


Interchange dispensers on the same system

- Non-contact picoliter to nanoliter to microliter
 - » 10 - 500 pL per drop
 - » 500 pL - 8 nL per drop
 - » 8 nL to 5 uL per dispense - single & multichannel
- Contact pin dispensing

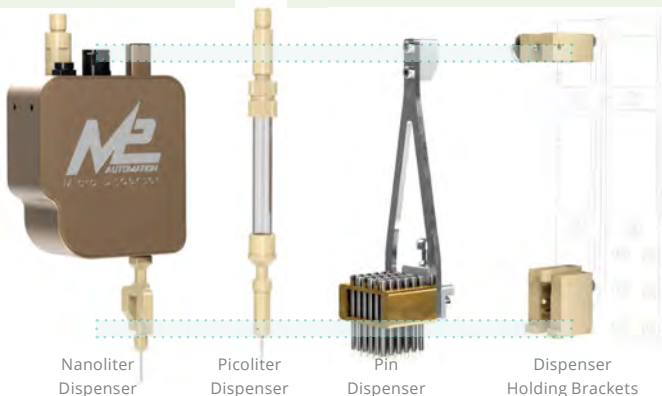

Common applications – supported by vision guidance

- Filling of microcavities, microwells
- Patterning lateral flow strips
- Coating biochip sensor pads
- Loading reagents in microfluidics devices
- Dispensing cells
- Spotting for multiplex immunoassays
- Printing protein, aptamer and DNA microarrays
- And so much more...

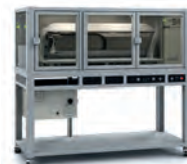


MODEL	iFOUR	iTWO	iONE - 600 / 1200	iZERO
Level	Compact Entry Level	R&D and Academic	Batch Production	In-Line Automated
Dimensions (W/D/H) & Weight	50 x 50 x 50 cm appr. 50 kg	110 x 75 x 110 cm appr. 150 kg	150/200 x 80 x 180 cm appr. 250 kg	85 x 151 x 248 cm appr. 500 kg
Built-in Windows PC	The built-in Windows PC with InDot Software is integrated as standard in all devices.			
Head-Camera	The Head Camera for in-line QC is integrated as standard in all devices.			
Washing Station	The Washing Station is integrated as standard in all devices.			
Travel Range	X 205 mm Y 200 mm Z 40 mm	X 250 mm Y 300 mm Z 100 mm	X 600 or 1200 mm Y 320 mm Z 44 mm	X 370 mm Y 398 mm Z 80 mm
Spottable Area	X 115 to 165 mm Y 195 mm	X 155 to 185 mm Y 300 mm	X 500 or 1100 mm Y 320 mm	X 370 mm Y 310 mm
Axis Resolution	< 0,5 µm	1 µm	1 µm	1 µm
Axis type	Belt-drives (XYZ)	Spindle drives (XYZ)	Linear Drives (XYZ)	Linear Drives (XYZ)
Position Repeatability	< +/-10 µm	+/-5 µm	+/-3 µm	+/-3 µm
Droplet Imaging	2D	2D	2D / 3D	2D / 3D
Humidity Control	✓	✓	✓	✓
Vacuum Table	✓	✓	✓	
Drying Station	✓	✓	✓	✓
HEPA Filter System with Flow Control		✓	✓	✓
De-Ionizer		✓	✓	✓
Cooling / Heating of source and target		✓	✓	✓
Orbital Shaker for source plates		✓	✓	✓
Ultrasonic Wash Station			✓	✓
On-the-fly dispensing			✓	✓
Analytical Balance for volume control			✓	✓
Conveyor / reel-to-reel				✓
Source Options 96, 384, 1536- MTP	1	2	Up to 10	Up to 4
3x8 postions for PCR tubes (100 or 200 µL)	1	2	Up to 10	Up to 4
Reagent vials of 0.5 to 1.5 mL	Up to 4	Up to 6	Up to 30	Up to 12
Custom holders	✓	✓	✓	✓
Target Options MTP Holder	2 plates	4 plates	10 or 22 plates	9 plate trays
Glass slide adapter for MTP holder (4 slides per MTP)	✓	✓	✓	✓
Planar Target Holder / Vacuum Table	Up to 10 slides (Depending on costumized target size)	Up to 28 slides (Depending on costumized target size)	Up to 176 slides (Depending on costumized target size)	Up to 65 slides (Depending on costumized target size)

iFOUR

The compact benchtop
Microdispenser



The entry level instrument for R&D applications

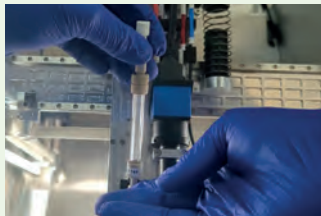

Easy-to-swap Dispenser
Options




Multiple source and target
options

iTWO

Designed for R&D and small
batch production



Optional high-duty mobile bench

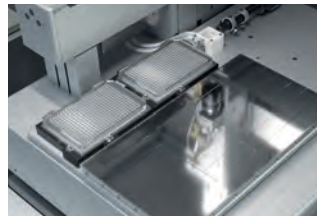

Cooling / Heating of source
and target



Various flat substrate
Vacuum Tables

iONE

The instrument for
high-production capacities



Available in two different sizes (iONE 600 or 1200)



High speed XYZ axes,
top speed: 1000 mm/s



Optional tray system for
ease-of-use

iZERO

The in-line-production
microarray-spotter system



Expandable modular production line to meet all needs




Integrated via Reel-to-Reel or
conveyor belts



Each module performs
partial spotting


Dispensers interchangeable on all Instruments

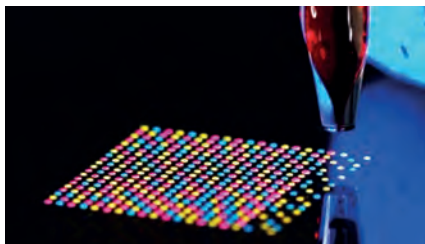


PDMD

PIEZO-DRIVEN

10 to 500 pL **NEW: 0.5 to 8 nL**

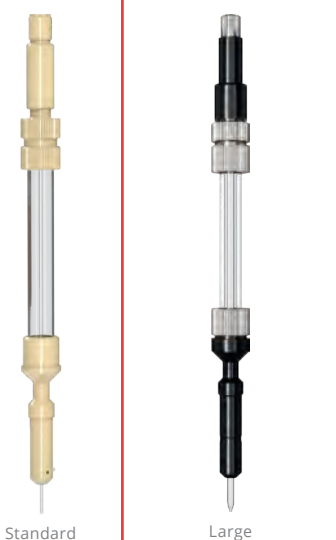


- ◆ 1 - 1000 Hz dispensing
- ◆ Volume CV < 2 %
- ◆ Different orifice diameters
- ◆ Inert borosilicate glass
- ◆ Robust construction for easy operator handling
- ◆ **NEW:** Up to 40 mg/mL protein dispensing capability

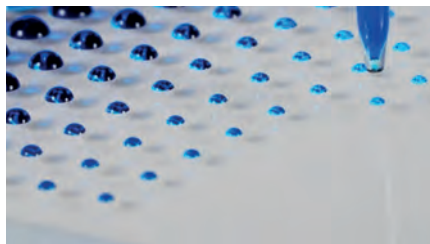


M2MD

VALVE-BASED

8 nL to 5 µL **NEW: 3 ch**

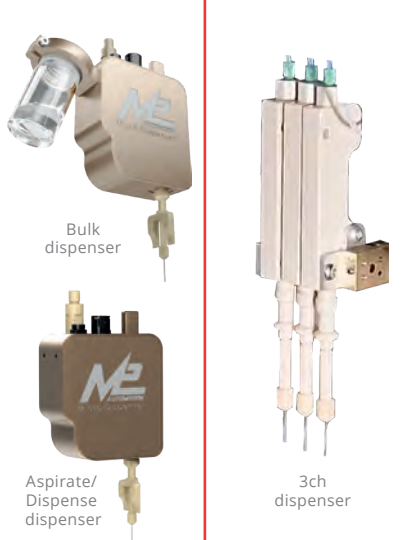

- ◆ Up to 50 Hz dispensing
- ◆ Volume CV 2 % to 5 %
- ◆ Bulk or Aspirate / Dispense
- ◆ Interchangeable tips:
 - » Glass tips
 - » Low volume plastic tips
 - » EP10 tips
- ◆ **NEW:** Multichannel standard & custom heads

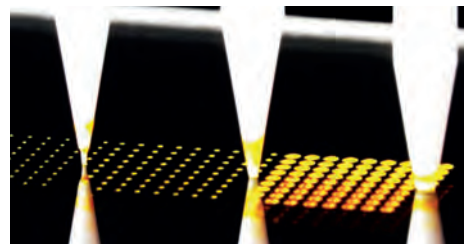

Pin-DMD

PIN-DRIVEN

350 pL to 35 nL **NEW: 1x1**



- ◆ Multiple configurations:
 - » 6x6, 2x8, 4x12, 1x1 ...
- ◆ Capillary or Solid (blunt)
- ◆ 80 µm to 500 µm tips
- ◆ Chemical inert ceramic tips
- ◆ Low cost metal pins available
- ◆ Matching wash stations
- ◆ **NEW:** single pin for quick tests



M2-Automation GmbH
Bessemerstraße 16
12103 Berlin
Germany



m2-automation.com
info@m2-automation.com
+49 30 856 11 939 0

For North and South America sales, service and support:
axiVEND LLC
350 E Crown Point Rd, Suite 1020
Winter Garden, FL 34787
United States



axivend.com
info@axivend.com
+1 833 294 8363

